


## general specification

The Amico Regal Series Surface Mounted Flatwall shall be manufactured by Amico Corporation, in accordance with job specific shop drawings and documents. The complete Flatwall assembly shall be cULus listed. The following is a general specification, and components listed may not be present or required in final product.

Figure 1. Regal Series 19" [483mm] Surface Mounted Flatwall

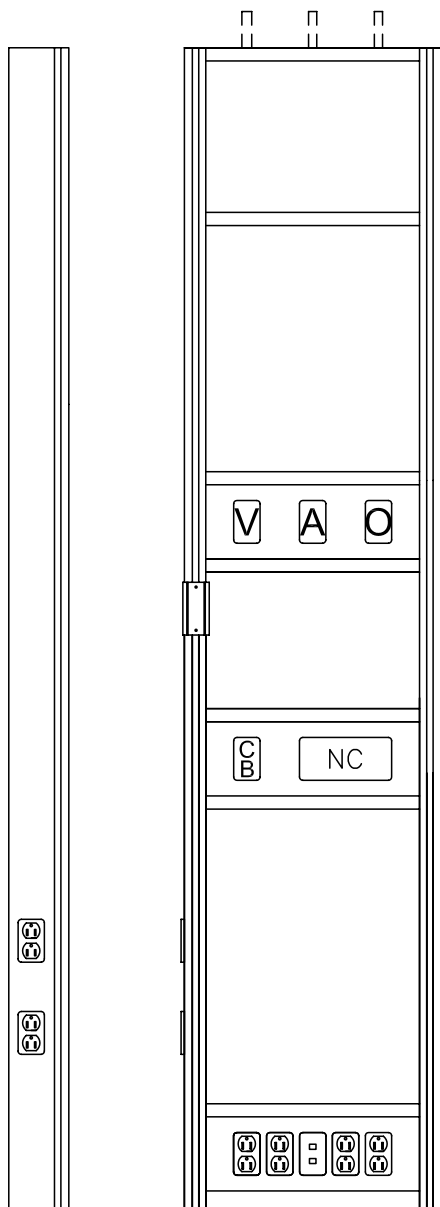

The Regal Series Flatwall System is listed by Underwriters Laboratories Inc.(UL) in Canada and the United States of America. This product also complies with the seismic requirements of, and is approved by the State of California's Office of Statewide Health Planning and Development (OSPHD). Amico Corporation manufactures all products in accordance with the National Fire Protection Association (NFPA), NFPA-99c (current edition) and National Electric Code (NEC).



project

---



---



---

The installation contractor is responsible for compliance with all local, state, and federal codes applicable to the installation of medical gas and electrical systems.

### submittals and verifications

Job specific shop drawings shall be produced for each project. These shop drawings will clearly indicate the area of Medical Gas termination, and electrical connection points inside or outside the Flatwall System. As Built drawings will be produced, and emailed at the time of shipment.

Amico Corporation will manufacture equipment as per signed approval drawings and verifications documents, which are provided in the job specific submittal package.

### materials and construction

Enclosure shall be constructed of extruded, anodized aluminum alloy sections to provide a modular surface mounted unit with integrated accessory rails for equipment management.

Standard sizes of the Regal Series Surface Mounted Flatwall are 19" [483mm], 28" [711mm], and 38" [965mm] widths and can be made to any custom height.

Each headwall system shall contain electrical conduit to enclose electrical wiring for each type of power (critical, normal & low voltage/communication). Electrical termination shall be located as indicated in the drawing.

All med gas piping shall be hard-piped and brazed to a single point of connection. Med gas pipes will terminate above the ceiling line.

All Fascia panels (aluminum and laminated) shall be removable for access for easy installation and maintenance of headwall services.

#### Integrated Accessory Rails:

The Regal Series Surface Mounted Flatwall has four (4) accessory channels, which are integrated into the aluminum extrusion assembly. The accessory channels are clear etched anodized.

#### Fascia:

Standard anodized aluminum finish where devices are located. Medium-Density-Fiberboard with High pressure laminate (HPL) elsewhere. Color and pattern to be selected by customer.

#### Medical Gas Outlets:

Amico Corporation manufactures several compatible connection types for Medical Gas and Vacuum Outlets. The outlet connection type shall be as called for on the signed approval documents.

#### Medical Gas Manifold:

Medical gas distribution shall be supplied by pipe drops to the single point connection above the unit as indicated on the drawing. All Medical Gas Outlets and piping shall be brazed and tested in accordance with NFPA 99c. All Medical Gas piping shall be Type L copper pipe. Each outlet, piping and manifold shall pass a 24 hour standing pressure test.

#### Electrical Wiring:

Wire for standard and critical branch power circuits shall be #10 or #12 (as specified) type THHN stranded

copper wire, 600 volt, with heat resistant thermo-plastic insulation for hot (black) and neutral (white). Grounds shall be #10 type R THHN stranded copper wire (green). All ground conductors shall be installed in conduit.

### **Grounding and Bonding:**

All ground conductors shall be installed in conduit. Each power receptacle shall have a ground conductor connected to a grounding screw. Where electrical terminations are located inside the flatwall, a grounding bus will be installed for each type of power, and shall be installed as to insure grounding for the complete power system.

### **Electrical Receptacles:**

Shall be Hubble Hospital Grade 15 or 20 amp, 120 or 277 volt, U.L listed and marked Hospital Grade. Quantity and type as shown on drawings. Amico shall furnish and install receptacles with wiring required, as shown on submittal shop drawings. All receptacles shall be Hubble hospital grade receptacles. Electrical contractor shall insure compatibility of plug on accessory equipment to be used with these devices:

### **Duplex receptacles shall be NEMA style 5-15R**

or 5-20R, color Ivory for use on normal power circuits, and color Red for use on critical branch power circuits.

**Simplex receptacles** shall be NEMA style 5-15R or 5-20R, color Ivory for use on normal power circuits, and color Red for use on critical branch power circuits.

**Safety receptacles** (if required) shall be duplex type, be NEMA style 5-15R or 5-20R, color Ivory for use on normal power circuits, and color Red for use on critical branch power circuits. Receptacles shall limit proper access to energized contracts and shall accept both 2 wire and 3 wire plugs.

**Locking receptacles** (if required) shall be simplex type, 20 amp, 120 or 277 volt, color black, and be Hubblelock 23000 HG style.

### **Low Voltage Data Provisions:**

Shall be connected to device junction box via conduit. Amico Corporation shall include pull cord extending from service provision to electrical termination point. Telephone jack, Blank faceplates, Data port or other Low Voltage device shall be furnished, installed and wired by Low Voltage contractor.

### **Switching:**

Shall be Hubble Industrial Grade 120 or 277 volt, 15 or 20 amps. Switch type options include SPST, 3-Way or Momentary, as shown on submittal shop drawings. Low Voltage Switching will be 0-12 volts, 15 amps unless otherwise noted. Amico Corporation shall furnish, pre-install and wire all switches.

### **Lighting (optional):**

Direct and Indirect lighting can be provided by the Amico 3' or 4' Overbed lighting system. Two fluorescent fixtures and ballasts shall be located in the Overbed light enclosure. These fixtures provide lighting on and above the patient bed. Each fixture shall accommodate four (4) fluorescent bulbs (supplied by others).

## **installation**

Mounting and final connections are the responsibility of the installing contractor. This includes receiving, storage, erection, overhead bracing, clean up, touch up, and carton disposal. The contractor shall supply all necessary installation materials, including such items as tools, fasteners, calking and electrical lamps.

The electrical contractor shall be responsible for the final electrical hook-up at service connection locations, as well as interconnection wiring on multi-sectional units. Once the installation is complete, the electrical contractor shall test equipment function and high voltage services / grounding in accordance with local electrical code requirements. It is the responsibility of the hospital and/or the approval authority to notify Amico if Local electrical code varies or differs from that of NEC. Amico manufactures all products in accordance with NEC (current edition).

The medical gas contractor shall be responsible for piping and final connection of all medical gas services. Amico Corporation will manifold medical gas plumbing internal to the Headwall System, and provide single point of connection as per the submittal documents. The medical gas contractor shall be responsible for purging, pressure testing, gas identification and system certification in accordance with NFPA 99c (current edition) or CSA-Z7396.01-06.

Please refer to the Regal Series Surface Mounted Flat-wall installation instructions for further detail.

## optional equipment

### **AMICO Accessories:**

Shall be provided in the quantities and types as ordered. Accessories shall be installed by others, in the locations determined by the end user. Please see the Amico Accessories catalog for additional product and mounting information.

## product code definition

### **FWXXXX-SM-CAR-P**

*12 Characters*

XXXX = Variable Width/Height (eg.1907)

CAR= Care Level (GEN, ICU, etc)

P = Private or Semi-Private (P, S)