

specifications

The medical gas outlet shall be an Amico Alert-1 series.

Outlets shall be manufactured with a 6 - 1/2" [165 mm] length type "K" 1/2" [12.7 mm] outside diameter (3/8" nominal) size copper inlet pipe stub, which is silver brazed to the outlet body. The body shall be of 1 - 5/16" [33 mm] diameter, one piece brass construction. For positive pressure gas services, the outlet shall be equipped with a primary and secondary check valve and the secondary check valve shall be rated at a maximum of 200 psi [1,379 kpa] in the event the primary check valve is removed for maintenance. Outlet bodies shall be gas specific by indexing each gas service to a gas specific dual pin indexing arrangement on the respective identification module.

A large color coded front plate shall be used for ease of gas identification and aesthetic appeal.

The latch valve assembly shall accommodate up to 3/4" [19 mm] wall thickness.


The latch-valve assembly shall be Ohmeda Diamond™ (Ohio Medical) Quick Disconnect compatible and accept only corresponding type gas specific adapters.

All M.R.I. outlets are manufactured from nonferrous materials.

All outlets shall be cleaned and degreased for medical gas service, factory assembled and tested.

The outlets are UL listed and CSA approved.

Amico products comply with NFPA-99.

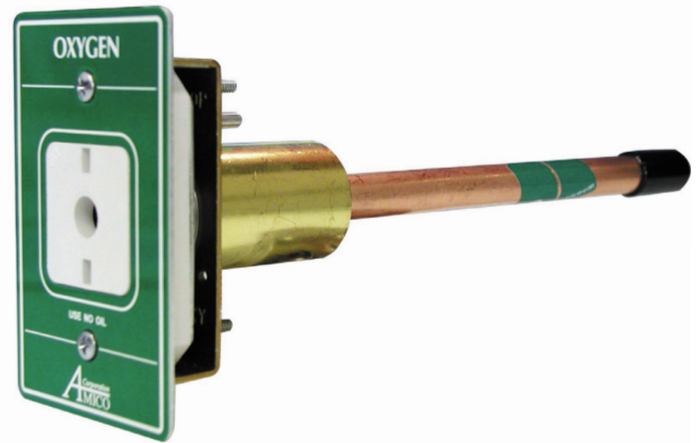


features

- ▶ Flow rate exceeds the requirements of NFPA and CSA
- ▶ Accepts only Ohmeda Diamond™ gas specific adapters
- ▶ 100 % Hydrostatic testing of each outlet and cleaned for medical gas service
- ▶ Indexed to eliminate interchangeability of gas services
- ▶ Up to 3/4" [19 mm] mounting plate thickness adjustment
- ▶ Universal Rough-In Assembly to accept Quick Disconnects or DISS Front Adapters (Can be interchanged at any time)
- ▶ Outlets are also available for M.R.I. use

project

ceiling column outlet quick disconnect - ohmeda compatible

Inch
[mm]

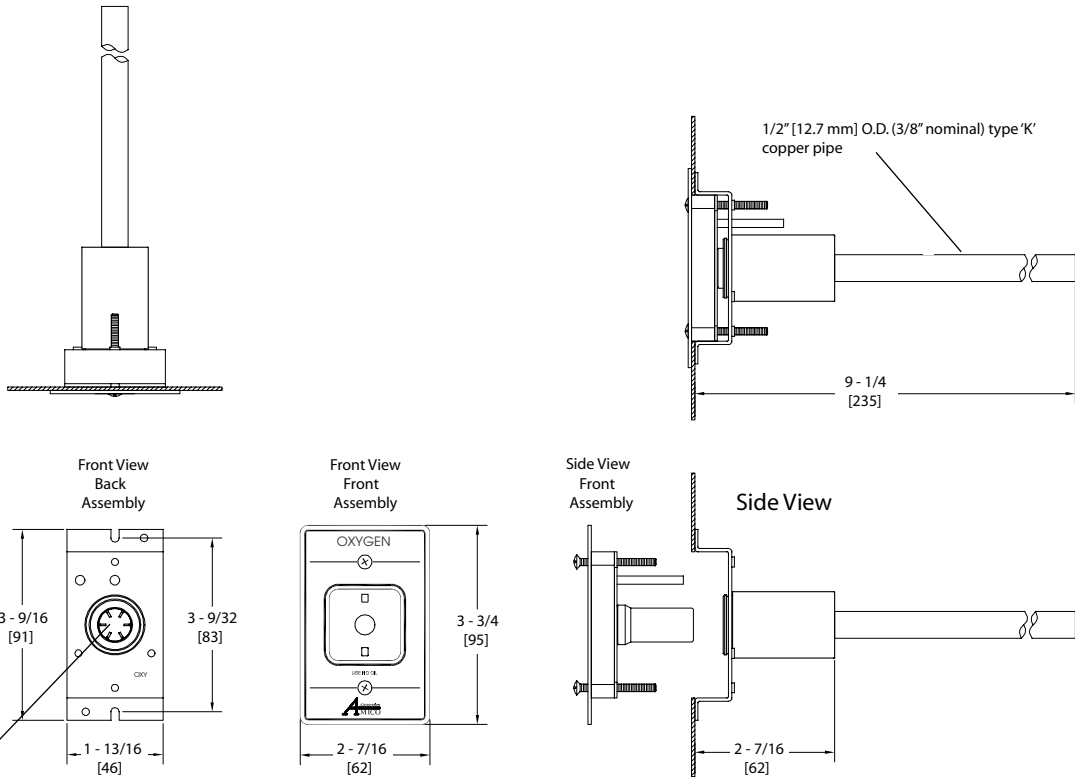

¹**NOTE:** For CSA/ISO, suction inlets are not supplied with a secondary check valve. A pressure plug is provided for testing purposes, rated at 150 psi [1,034 kPa].

model numbers

O-QDCCOL-L-XXX

The L defines the Language: _____

U = English (NFPA) Oxy, Air, Vac and WAG only
 E = English (CSA/ISO)
 F = French (CSA/ISO)
 S = Spanish (NFPA/ISO)

Example:

English NFPA, Oxygen = O-QDCCOL-U-OXY

The XXX defines the Gas: _____

Oxygen	= OXY
Medical Air	= AIR
Nitrous Oxide	= N2O
MedVac ¹	= VAC
Waste Anesthetic Gas Disposal (NFPA)	= WAG
Anesthetic Gas Scavenging System ¹ (CSA)	= AGS
Carbon Dioxide	= CO2

¹Please see note above

represented by: