

SPECIFICATIONS:

The 37-pin connector shall be an Amico Alert-2 series.

The kit provides a dedicated 37-pin bed port for communication between beds and controlled devices (i.e. Nurse calls, LVC, TV and Radio).

Pin No	Color Code	Typical Pin Functions
Cable "A"		
2	Blk	Direct Light, Right
3	Wht	Indirect Light, Right
20	Red	Direct Light, Left
21	Grn	Indirect Light, Left
27	Brn	Light Common
Cable "B"		
4	Yel	Speaker High, Right
7	Orn	Nurse Call Interlock
8	Blu	Audio Transfer (-)
9	Brn	Audio Transfer (+)
22	Grn	Speaker High, Left
25	Red	Nurse Call (+)
26	Wht	Nurse Call (-)
35	Blk	Speaker Low Common
Cable "C"		
5	Blu/Blk	Volume Ctrl. Wiper, Right
6	Red/Blk	Radio Control, Right
10	Blk/Red	Interlock (+)
11	Wht/Red	Interlock (-)
14	Blk/Wht	Volume Ctrl. Low Common
15	Red/Wht	Volume Ctrl. High Common
16	Orn	Red Nurse Ans. Light (+)
18	Red	D.S. (+)
19	Grn	Yellow Nurse Call Light (+)
23	Grn/Wht	Volume Control. Left
24	Orn/Blk	Radio Control, Left
28	Blk	Yellow Nurse Call Light (-)
29	Wht	Red Nurse Ans. Light (-)
30	Orn/Red	Bed Exit Priority Call (+)
31	Blu/Wht	Bed Exit Priority Call (-)
33	Wht/Blk	TV Channel Stepping (-)
34	Blu	TV Channel Stepping (+)
37	Grn/Blk	Radio Control, Common

Exact pin use depends on function and type of bed being used. Not all beds use every pin shown.




FEATURES:

- Stainless Steel or Aluminum wall plate with 37-pin connector and Plug.
- 12 foot cables for connections in headwalls or bed locators.
- Each cable is stripped back 3" and then each conductor is stripped 1".
- Custom lengths available on special order.

MODEL NUMBERS

W-37PIN-SS-12

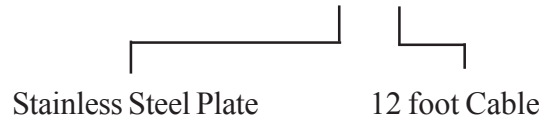

PROJECT:

D I S T R I B U T E D B Y :

A *Corporation*
AMICO 14 Madison Rd., Fairfield, NJ 07004 USA
85 Fulton Way, Richmond Hill, ON L4B 2N4 Canada

Tel: 1-877-GOAMICO (1-877-462-6426)
Fax: (905) 764-0862 www.amico.com

JUL-04